

Protein A HP SpinTrap Protein A HP MultiTrap Protein A/G SpinTrap Buffer Kit

Protein A HP SpinTrap™ and Protein A HP MultiTrap™ (Fig 1) are designed for the enrichment of proteins of interest, immunoprecipitation, and the purification of monoclonal and polyclonal antibodies from cell lysates and biological fluids. Protein A HP SpinTrap are prepacked, single-use spin columns. Protein A HP MultiTrap are prepacked 96-well filter plates. Both products contain Protein A Sepharose™ High Performance.

The spin columns and 96-well filter plates address the need for flexible, small-scale preparation of protein samples before downstream protein analyses such as gel electrophoresis, liquid chromatography, and LC-MS. When using the columns and plates for protein enrichment, an antibody of choice is first immobilized onto the Protein A media, and then used for the capture and enrichment of the protein of interest. The protein can be enriched several hundredfold, depending on the specificity of the antibody.

The benefits of Protein A HP SpinTrap and Protein A HP MultiTrap when running the protein enrichment/immunoprecipitation protocols are:

- Reproducible performance from run to run enables quantitative and comparative expression studies
- High performance for yield and purity, which can be further increased using the supplied Optimization Guide that supports performance optimization, protocol modifications, and troubleshooting.
- Protocols that have been tested and optimized for downstream analyses such as electrophoresis and LC-MS.



Fig 1. Protein A HP SpinTrap columns and MultiTrap 96-well filter plates are designed for efficient, small-scale enrichment of proteins/immunoprecipitation from cell lysates and biological fluids and can be used with the Protein A/G HP SpinTrap Buffer Kit for increased convenience. They can also be used for antibody purification.

- Choice of protocols; a classic protocol for speed or a cross-link protocol ensuring elution of the protein of interest only, with the antibody remaining bound to the column.

For antibody purification, the benefits of Protein A HP SpinTrap and Protein A HP MultiTrap are:

- High purity and yield of antibodies
- Simple and proven methods for reproducible results
- Quick purification due to the prepacked format and the absence of pretreatment.



Buffer Kit

Protein A/G SpinTrap Buffer Kit (Fig 1) is designed for protein enrichment using Protein A HP or Protein G HP SpinTrap columns. The kit is sufficient for 16 reactions and contains reagents for both the Classic and Cross-link protocols (with or without elution of antibody together with target protein). Some of the reagents are delivered in two bottles for increased stability. The kit consists of both stock solutions and reagents ready for use. Working solutions are prepared by adding distilled water directly into the stock solution bottle.

The kit eliminates time-consuming buffer preparation and thus promotes fast, reproducible and convenient enrichment of target protein from a complex protein sample.

Characteristics

Protein A HP SpinTrap and Protein A HP MultiTrap are both prepacked with Protein A Sepharose High Performance, a proven medium with strong affinity for IgG subclasses. Reliable sample preparation is achieved through fast kinetics and the high binding capacity of the medium. Specificity for capture of the protein is achieved by immobilization of antibodies to the Sepharose beads through well-established coupling techniques (Fig 4). The agarose-based medium provides a hydrophilic and chemically favorable environment for coupling, while the highly cross-linked structure of the 34- μ m spherical beads ensures excellent flow of sample through the spin columns and the 96-well filter plates. Table 1 lists the main characteristics of the products.

Table 1. Characteristics of prepacked Protein A HP SpinTrap columns, Protein A HP MultiTrap 96-well filter plates and the Protein A/G Buffer Kit.

Prepacked medium	Protein A Sepharose High Performance
Matrix	Highly cross-linked agarose, 6%
Ligand	Native protein A
Ligand coupling method	N-hydroxysuccinimide activation
Ligand density	Approx. 3 mg protein A/ml medium
Binding capacity ¹	> 10 mg human IgG/ml medium
Average particle size	34 μ m
pH stability ²	3–9 (long term), 2–9 (short term)
Working temperature	4°C to 30°C
Storage solution	20% ethanol
Storage temperature	4°C to 30°C

Protein A HP SpinTrap

Volume of prepacked medium	100 μ l
Column volume	800 μ l
Column material	Polypropylene and polyethylene

Protein A HP MultiTrap

Filter plate size ³	127.8 × 85.5 × 30.6 mm
Prepacked medium volume/well	50 μ l
Well volume	800 μ l
Filter plate material	Polypropylene and polyethylene
Centrifugation speed ⁴	700 × g
Vacuum pressure	
Recommended	-0.1 to -0.3 bar
Maximum	-0.5 bar

¹ Protein dependent

² pH below 3 is sometimes required to elute strongly bound antibody species. However, protein ligands may hydrolyze at very low pH

³ According to American National Standards Institute (ANSI) and Society for Biomolecular Screening (SBS) standards 1-2004, 3-2004, and 4-2004

⁴ Actual settings depend on the sample properties and pretreatment

Protein A/G SpinTrap Buffer Kit

Buffer	Content	Formulation	Volume
Binding/Washing Buffer (TBS)	0.5 M Tris; 1.5 M NaCl, pH 7.5	10x	2x5 ml
Elution Buffer (Classic)	2.5 % HAc	Ready to use	20 ml
Elution Buffer (Cross-link)	1 M Glycine-HCl, pH 2.9	10x	3 ml
Cross-link Solution A	2 M Triethanolamine, pH 8.9	10x	4 ml
Cross-link Solution B	1 M Ethanolamine, pH 8.9	10x	2x1 ml

Antibody purification

The purification of antibodies with Protein A HP SpinTrap and MultiTrap can be divided into four stages; equilibration, sample incubation, wash, and elution. Each step involves centrifugation.

Binding of the antibody is performed at neutral pH, and elution by lowering the pH. Eluted material is collected in tubes containing neutralizing buffer to preserve the activity of acid-labile IgG's.

Purifying antibodies with Protein A HP SpinTrap or MultiTrap is a simple, four-stage procedure that can be completed in less than 20 min: (1) Equilibrate by adding binding buffer and centrifuge; (2) add sample and incubate; (3) wash with binding buffer; (4) elute the target protein with elution buffer.

Recommended buffers:

Binding buffer: 20 mM sodium phosphate, pH 7.0

Elution buffer: 0.1 M glycine-HCl, pH 2.7

Neutralizing buffer: 1 M Tris-HCl, pH 9.0

Ab Buffer Kit is available as an accessory for increased convenience. The kit contains buffers for Binding, Elution and Neutralizing, optimized for rapid purification of monoclonal and polyclonal IgG with immobilized Protein A or Protein G media. The kit eliminates time-consuming buffer preparation and thus promotes fast, reproducible and convenient purification work.

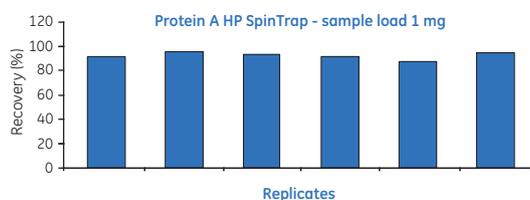


Fig 2. Antibody purification of human IgG (Gammanorm), 1 mg loaded/ column, performed according to Ab purification protocol.

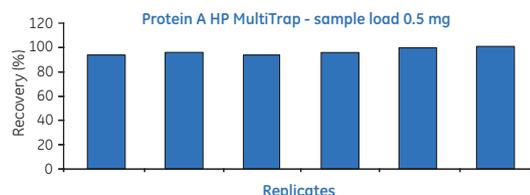


Fig 3. Antibody purification of human IgG (Gammanorm), 0.5 mg loaded/ multiplate well, performed according to Ab purification protocol using buffer volumes adjusted for 50 µl media/well.

Table 2. Relative binding strengths of antibodies from various species to protein A as measured in a competitive ELISA test. Antibody binding strength to protein G is also shown for comparison. GE Healthcare also offers prepacked Protein G HP SpinTrap and Protein G HP MultiTrap; see Ordering information for details.

Species	Subclass	Protein A binding	Protein G binding
Human	IgA	variable	—
	IgD	—	—
	IgE	—	—
	IgG ₁	++++	++++
	IgG ₂	++++	++++
	IgG ₃	—	++++
	IgG ₄	++++	++++
	IgM*	variable	—
Avian egg yolk	IgY†	—	—
Cow		++	++++
Dog		++	+
Goat		—	++
Guinea pig	IgG ₁	++++	++
	IgG ₂	++++	++
Hamster		+	++
Horse		++	++++
Koala		—	+
Llama		—	+
Monkey (rhesus)		++++	++++
Mouse	IgG ₁	+	++++
	IgG _{2a}	++++	++++
	IgG _{2b}	+++	+++
	IgG ₃	++	+++
	IgM*	variable	—
Pig		+++	+++
Rabbit		++++	+++
Rat	IgG ₁	—	+
	IgG _{2a}	—	++++
	IgG _{2b}	—	++
	IgG ₃	+	++
Sheep		+/-	++

* Purified using HiTrap™ IgM Purification HP columns, code no. 17-5110-01

† Purified using HiTrap IgY Purification HP columns, code no. 17-5111-01

++++ = strong binding

++ = medium binding

— = weak or no binding

Reproducible protein enrichment/immunoprecipitation

Enrichment of a particular protein is often desired to increase its signal in subsequent analysis steps. Protein A HP SpinTrap and MultiTrap are used to immobilize a biospecific ligand (i.e. an antibody) with affinity for the protein of interest. The Protein A/G SpinTrap Buffer Kit is optimized for the enrichment protocol and improves reproducibility and convenience.

To correlate protein expression with, for example, disease or treatment, large numbers of samples must be prepared and analyzed. This is tedious work and a source of error. SpinTrap columns and MultiTrap 96-well filter plates provide reproducibility, flexibility, and convenience, thus minimizing variation during the enrichment step.

One of the major advantages of Trap products are the protocols which contain full descriptions of the components of each product. They also contain an optimization guide and recommended elution buffers to enable the direct analysis of the eluates with electrophoresis or LC-MS. SpinTrap columns require only a standard microcentrifuge. MultiTrap 96-well filter plates allow sample preparation by centrifugation or vacuum, either operated manually or automated using robotics.

Specific coupling of antibodies

Protein A Sepharose High Performance has high affinity for the Fc region of antibodies in a variety of species (Table 2). Note that Protein G has affinity for a wider range of antibody species than Protein A. The protein enrichment protocols provided offer two methods for antibody attachment, allowing the enriched protein to be either eluted together with the antibody (classic protocol) or without (cross-link protocol). Elution can be performed using the buffers described in Table 3.

Table 3. Suggested elution buffers for various situations in protein enrichment/immunoprecipitation ¹

Glycine/HCl, 1-2 M urea, pH 2.5-3.1	Most antibody-antigen bonds are broken and sufficient elution is often achieved. This is often a first-choice buffer to screen for optimal elution conditions.
Glycine/HCl, pH 2.5-3.1	Many antibody-antigen bonds are broken and sufficient elution is often achieved.
0.5 M acetic acid	Low pH buffer compatible with mass spectrometry due to the volatility of acetic acid.
2% SDS	Breaks all protein-protein bonds and solubilizes even the most difficult proteins. Can be used in aqueous solution or as an additive to other buffers. SDS is often a constituent of electrophoresis loading buffer making it compatible with many electrophoresis procedures.
Citric acid, pH 2.5-3.1	Many antibody-antigen bonds are broken and sufficient elution is often achieved. This buffer has performed well when used in the classic protocol in our labs.
0.1 M ammonium hydroxide	A basic elution buffer used, for instance, when the protein of interest is acid labile.

¹ For further suggestions regarding protocol optimization, see Instructions 28-9067-70 (Protein A HP SpinTrap) and 28-9067-71 (Protein A HP MultiTrap)

Easy-to-use protocols enable reproducible protein enrichment

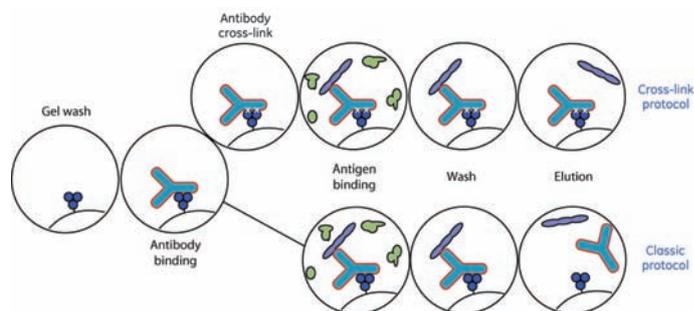


Fig 4. Schematic representation of the alternative cross-link and classic protocol approaches for use with Protein A HP SpinTrap and Protein A HP MultiTrap.

Cross-link protocol for protein enrichment

The major advantage of the cross-link protocol is that the antibody is covalently bound to the Protein A Sepharose High Performance medium, which enables the elution of the protein without co-elution of the antibody. This is often advantageous since the antibody is generally in excess compared to the protein of interest. High levels of antibody in the eluted fractions can obscure the desired signal from the protein of interest.

To demonstrate the functionality of the cross-link protocol using Protein A HP SpinTrap, human transferrin was enriched from a background of *E. coli* protein sample¹. The transferrin concentration was 0.15% of the total *E. coli* protein content, which approximately corresponds to the concentration of a medium-abundant protein. Capture of the protein of interest was achieved using polyclonal rabbit anti-human transferrin that was immobilized to the protein A ligand.


As shown in Figure 5, the majority of the enriched protein was eluted in the first elution step. In proteomics applications, it is often desirable to achieve the highest possible concentration of the protein of interest since sensitivity and detection limits are frequently limiting factors.

Reproducibility of Protein A HP SpinTrap is indicated by the analysis of the recovery (percentage of start material) of three replicates that were run in parallel (Fig 5B). Quantitation of eluted protein of interest was performed using a standard curve with known amounts of transferrin (data not shown). Essentially no unwanted protein was bound to the medium after five washing steps, as indicated by the fact that very low amounts of protein were detected in the fifth washing step (Fig 5A, lane 5). In this particular

¹ In all examples described in this Data File, model proteins were enriched after being added to *E. coli* protein sample.

setup, an enrichment in the range of 200-fold was achieved relative to the start material using the cross-link protocol. Variation between spin columns (relative standard deviation) is generally below 15% with respect to purity and below 10% with respect to recovery (see Data File 28-9067-91 AA, Streptavidin HP SpinTrap and Streptavidin HP MultiTrap).

Trap product: Protein A HP SpinTrap
Sample: 5 mg/ml *E. coli* protein containing 7.5 µg/ml human transferrin
Sample volume: 0.2 ml
Antibody: Polyclonal rabbit anti-human transferrin
Binding buffer: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5)
Wash buffer: TBS, 2 M urea, pH 7.5
Elution buffer: 0.1 M glycine/HCl, 2 M urea, pH 3

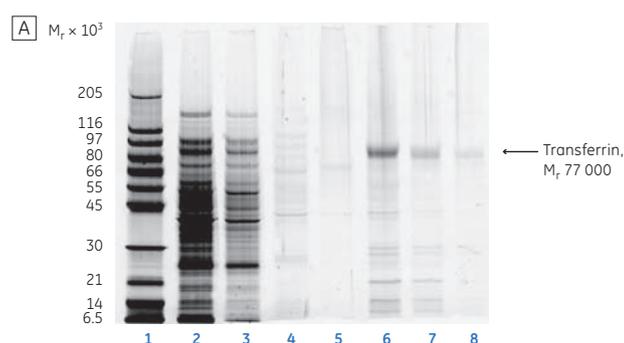

Lanes

1. Protein Molecular Weight Standards (broad range)
2. Flowthrough (diluted 1:30)
3. First wash (diluted 1:10)
4. Third wash
5. Fifth wash
6. First elution
7. Second elution
8. Third elution


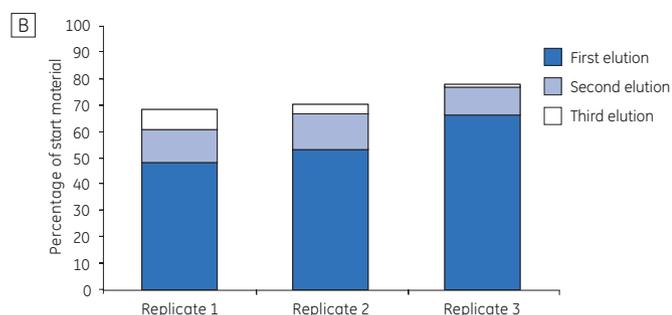


Fig 5. Enrichment of transferrin from *E. coli* cell lysate using the cross-link protocol for Protein A HP SpinTrap. Fractions were collected from every step of the enrichment process and analyzed by **(A)** SDS-PAGE (wash steps 2 and 4 have been omitted from the gel). The result from one of the three replicates is shown. The gel was poststained with Deep Purple Total Protein Stain and scanned using Ettan™ DIGE Imager. **(B)** Each elution step was quantitated using a standard curve of known amounts of transferrin. Quantitation of each elution step was performed in triplicate.

The cross-link protocol also enables highly reproducible results on Protein A HP MultiTrap. Analysis of 10 parallel replicates in the enrichment of human serum albumin (HSA) from *E. coli* lysate shows the reproducibility of the cross-link protocol used with Protein A HP MultiTrap (Fig 6). Well-to-well variation was below 10% (relative standard deviation) for both purity and recovery.

Trap product: Protein A HP MultiTrap
Sample: 5 mg/ml *E. coli* protein containing 7.5 µg/ml HSA
Sample volume: 0.2 ml
Antibody: Polyclonal rabbit anti-human albumin
Binding buffer: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5)
Wash buffer: TBS, 2 M urea, pH 7.5
Elution buffer: 0.1 M glycine/HCl, 2 M urea, pH 3

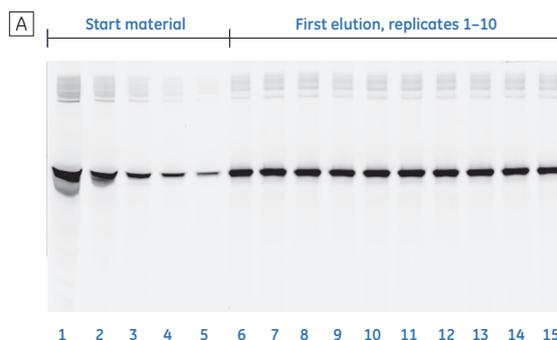

Lanes

1. Start material
2. Start material (diluted 1:2)
3. Start material (diluted 1:4)
4. Start material (diluted 1:8)
5. Start material (diluted 1:16)
- 6-15. First elution step, 10 replicates

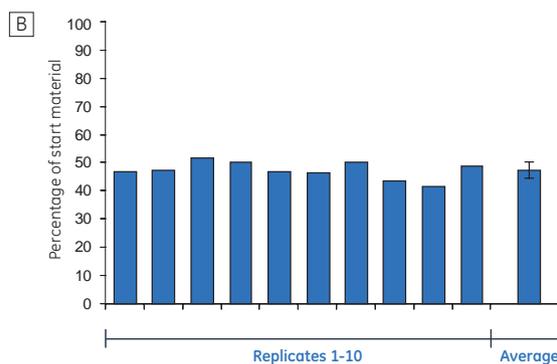



Fig 6. Enrichment of HSA from *E. coli* cell lysate using the cross-link protocol for Protein A HP MultiTrap. **(A)** Analysis by SDS-PAGE of 10 replicates of the first elution step. For visualization purposes, the HSA antigen was labeled with CyDye DIGE Fluor Cy5 minimal dye. After electrophoresis, the gel was scanned directly using Ettan DIGE Imager. **(B)** Percentage recovery relative to the start material in the first elution step. The average value is indicated to the right in the graph and the relative standard deviation of the 10 replicates is indicated by the error bar.

Classic protocol for protein enrichment

The major advantage of the classic protocol is that it is very convenient and fast since no coupling reaction needs to be performed. The price paid for the time-saving aspect of the classic protocol is co-elution of the antibody with the protein of interest. In many cases this is acceptable, for example, when sample is labeled with CyDye™ or radiolabeled. In these instances, only labeled proteins will be visualized through the respective detection procedures leaving IgG molecules undetected. Other instances when the classic protocol can be applied are when the protein of interest is of significantly different size than IgG, enabling the separate analysis of the protein of interest through, for example, electrophoresis.

To demonstrate the efficiency of the classic protocol, HSA was enriched from a background of *E. coli* protein sample (concentration of protein of interest was 0.15%) using Protein A HP SpinTrap.


Trap product: Protein A HP SpinTrap
 Sample: 5 mg/ml *E. coli* protein containing 7.5 µg/ml HSA
 Sample volume: 0.2 ml
 Antibody: Polyclonal rabbit anti-human albumin
 Binding and wash buffer: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5)
 Elution buffer: 0.5 M acetic acid

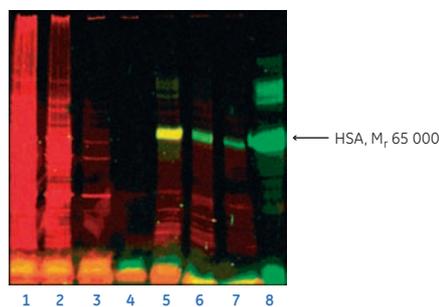

Lanes
 1. Flowthrough (diluted 1:30)
 2. First wash (diluted 1:10)
 3. Third wash
 4. Fifth wash
 5. First elution
 6. Second elution
 7. Third elution
 8. Pure HSA (1.5 µg, Cy3 labeled)

Fig 7. Enrichment of CyDye labeled (CyDye DIGE Fluor minimal dye) HSA from *E. coli* lysate using the classic protocol with Protein A HP SpinTrap. Fractions were collected from every step of the enrichment process and analyzed by SDS-PAGE (wash steps 2 and 4 have been omitted from the gel). For visualization purposes, *E. coli* protein was labeled with Cy3 minimal dye and HSA was labeled with both Cy3 and Cy5 minimal dyes. After electrophoresis, the gel was scanned directly using Ettan DIGE Imager. The antibody was not labeled and is therefore not visible in the image.

Trap product: Protein A HP MultiTrap
 Sample: 5 mg/ml *E. coli* protein containing 7.5 µg/ml human transferrin
 Sample volume: 0.2 ml
 Antibody: Polyclonal rabbit anti-human transferrin
 Binding and wash buffer: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5)
 Elution buffer: 0.5 M acetic acid

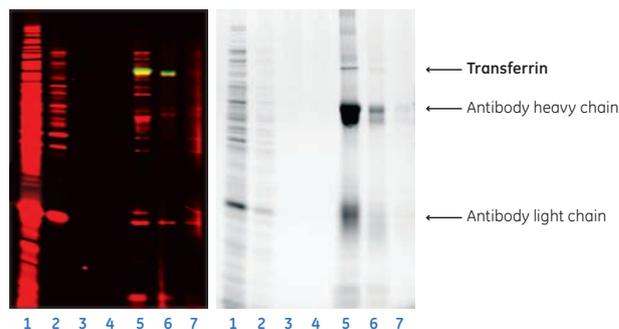

Lanes
 1. Flowthrough (diluted 1:30)
 2. First wash (diluted 1:10)
 3. Third wash
 4. Fifth wash
 5. First elution
 6. Second elution
 7. Third elution

Fig 8. Enrichment of transferrin from *E. coli* cell lysate using the classic protocol for Protein A HP MultiTrap. Fractions were collected from every step of the enrichment process and analyzed by SDS-PAGE (wash steps 2 and 4 have been omitted from the gel). For visualization purposes the *E. coli* protein was labeled with CyDye DIGE Fluor Cy5 minimal dye while the transferrin was labeled with both Cy3 and Cy5 minimal dyes. After electrophoresis, the gel was scanned directly using Ettan DIGE Imager (left image). The antibody was not labeled and is therefore not visible in this image. The gel was also poststained with Deep Purple Total Protein Stain to visualize total protein including antibodies (right image).

Analysis by SDS-PAGE of collected fractions from the spin column run revealed a significant enrichment of the protein (Fig 7).

Figure 8 shows analysis by SDS-PAGE of enrichment of transferrin using the classic protocol with Protein A HP MultiTrap. In addition to the Cy3™ and Cy5™ scanning, the gel was poststained with Deep Purple™ Total Protein Stain revealing the eluted antibodies on the gel. Transferrin is larger than IgG heavy chain (reducing conditions), which enabled detection by total protein-staining techniques.

Although the classic protocol is fast and convenient, it is more sensitive to harsh washing conditions compared to the cross-link protocol. This is due to the fact that the interaction between the Fc part of the antibody and protein A is sensitive to urea and acidic conditions. In this particular experiment, only Tris buffered saline (TBS), pH 7.5 was used as wash buffer.

Trap product: Protein A HP MultiTrap
 Sample: 5 mg/ml *E. coli* protein containing 7.5 µg/ml HSA
 Sample volume: 0.2 ml
 Antibody: Polyclonal rabbit anti-human albumin
 Binding and wash buffer: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5)
 Elution buffer: 0.5 M acetic acid

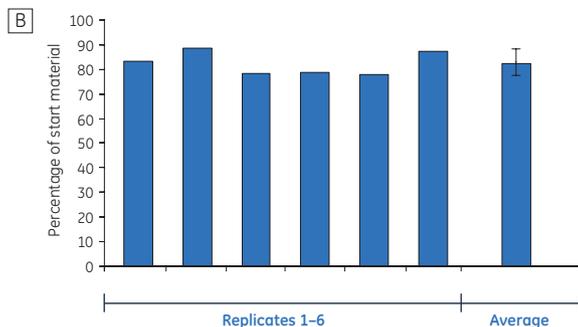
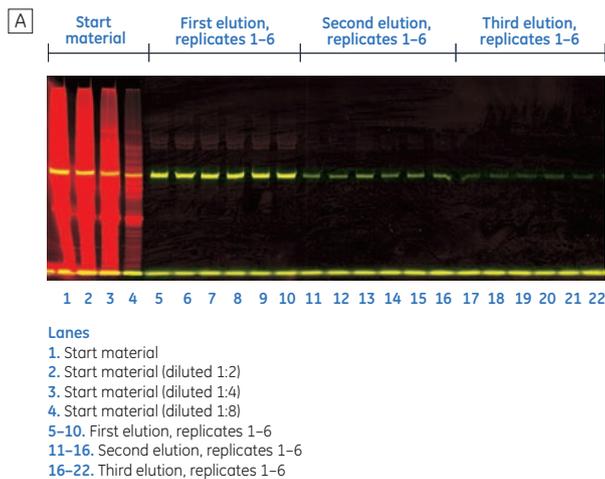


Fig 9. Enrichment of HSA from *E. coli* cell lysate using the classic protocol with Protein A HP MultiTrap. Six replicates of the enrichment are shown. Fractions were collected from the three elution steps and analyzed by **(A)** SDS-PAGE. For visualization purposes, the *E. coli* protein was labeled with CyDye DIGE Fluor Cy5 minimal dye and HSA was labeled with both Cy3 and Cy5 minimal dyes. After electrophoresis, the gel was scanned directly using Ettan DIGE Imager. The antibody was not labeled and is therefore not visible. **(B)** Average recovery of HSA after enrichment.

An additional advantage of the classic protocol compared to the cross-link protocol is that a larger part of the protein of interest is often found in the early elution fractions. This is exemplified in Figure 9 where approximately 80% of the total recovery is found in the first elution step. In concurrence with the cross-link protocol (Fig 4), a high level of reproducibility is obtained with Protein A HP MultiTrap; a relative standard deviation of less than 5% with respect to recovery in the first elution step was observed.

Label-independent enrichment in protein analysis workflows

In protein analysis workflows, CyDye labeling techniques are often used to enable protein quantitation and minimize experimental variation. An investigation was therefore performed to determine the usability of Protein A HP SpinTrap when comparing Cy3 and Cy5 differentially labeled protein samples.

In the study, human transferrin, labeled with CyDye DIGE Fluor Cy3 and Cy5 minimal dyes, was enriched from *E. coli* lysate. Labeled transferrin was added to the lysate at various ratios ranging from 1:3 to 2:1 with respect to Cy3 and Cy5 labels. The protein of interest was enriched according to the classic protocol and Cy3: Cy5 ratio differences were analyzed and quantitated (Fig 10). Expected levels of Cy3 and Cy5 labeled protein corresponded well to the measured levels of labeled protein ($R^2 > 0.99$) demonstrating that Protein A SpinTrap performed equally well for both the Cy3 and Cy5 labeled protein of interest (Fig 10B).

Trap product: Protein A HP SpinTrap
 Sample: 5 mg/ml *E. coli* protein containing various ratios of Cy3: Cy5 labeled human transferrin
 Sample volume: 0.2 ml
 Antibody: Polyclonal rabbit anti-human transferrin
 Binding buffer: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5)
 Wash buffer: TBS, pH 7.5
 Elution buffer: 0.1 M glycine/HCl, pH 2.8

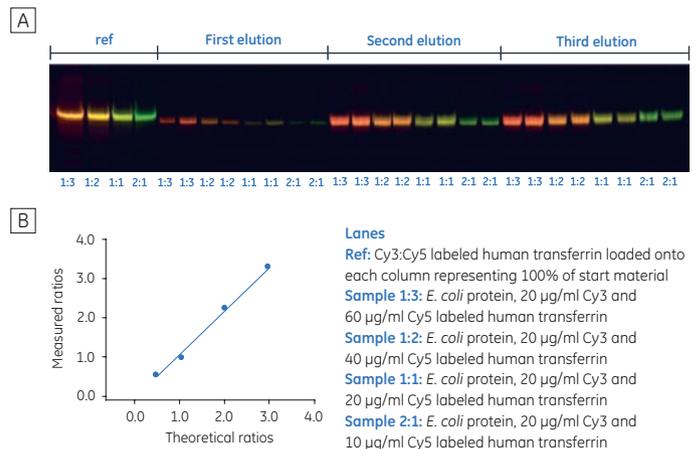


Fig 10. Enrichment of transferrin labeled with CyDye DIGE Fluor Cy3 or Cy5 minimal dye from *E. coli* cell lysate at different ratios ranging from 1:3 to 2:1 with respect to Cy3 and Cy5 labeled protein. Protein A HP SpinTrap was used applying the classic protocol for the enrichment. **(A)** Samples were collected from the three elution steps and analyzed by SDS-PAGE. After electrophoresis, the gel was scanned directly using Ettan DIGE Imager (antibody and *E. coli* protein were not labeled and are therefore not visible). **(B)** Quantitation of Cy3 and Cy5 labeled transferrin enriched on the SpinTrap columns was performed using ImageQuant™ TL software. Measured values in the first elution steps were compared to theoretical values.

Improved recovery of low-abundant protein

Products: Protein A HP SpinTrap and Pierce Seize Classic (A)
Sample: 5 mg/ml *E. coli* protein containing 7.5 µg/ml human transferrin
Sample volume: SpinTrap: 100 µl medium/column, Seize: 200 µl medium/column
Antibody: Polyclonal rabbit anti-human transferrin
Binding and wash buffer: SpinTrap: Tris buffered saline (TBS: 50 mM Tris, 150 mM NaCl, pH 7.5); Seize: BupH Modified Dulbecco's PBS pack
Elution buffer: SpinTrap: 2.5% HAC; Seize: Immunopure IgG Elution buffer
SDS-PAGE sample buffer: SpinTrap: 4* NSB; Seize: Non-reducing Lane Marker Sample Buf


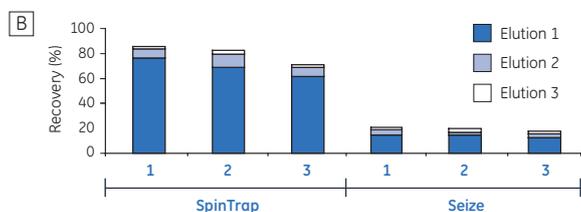


Fig 11. Comparison of enrichment of transferrin by Protein A HP SpinTrap and Pierce Seize Classic (A). The standard protocols and recommended buffers for the respective products were applied. Transferrin labeled with Cy5 was enriched from *E. coli* cell lysate. Protein A HP SpinTrap was used with the classic protocol for enrichment. **(A)** Samples were taken from the three elution steps and analyzed by SDS-PAGE. After electrophoresis, the gel was scanned directly using Ettan DIGE Imager (antibody was not labeled and is therefore not visible). **(B)** Approximately four times higher recovery was achieved by Protein A HP Spin Trap compared with Pierce Seize Classic (A). The purity at $37 \pm 4\%$ was similar for both products.

For contact information for your local office, please visit, www.gelifesciences.com/contact

www.gelifesciences.com/trap

GE Healthcare Bio-Sciences AB
 Björkgatan 30
 751 84 Uppsala
 Sweden

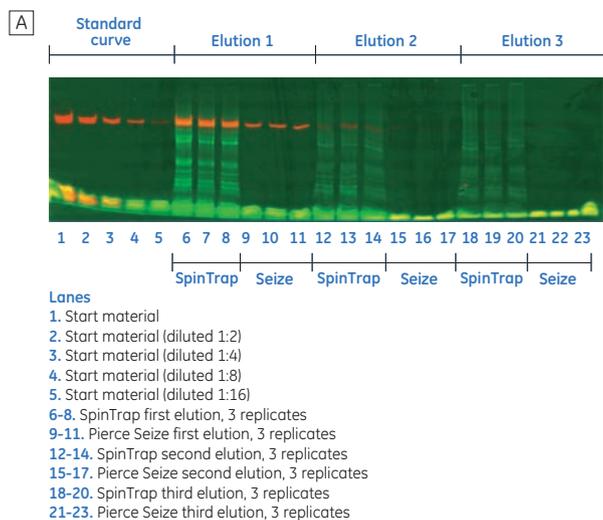

Ordering information

Products	Quantity	Code no.
Protein A HP SpinTrap	16 columns	28-9031-32
Protein A HP MultiTrap	4 × 96-well plates	28-9031-33
Protein A/G HP SpinTrap Buffer Kit	1	28-9135-67
Collection plate 500 µl V-bottom (for collection of fractions from MultiTrap)	5 × 96-well plates	28-4039-43
Related products		
Sample Grinding Kit	50 samples	80-6483-37
Protease Inhibitor Mix	1 ml	80-6501-23
Nuclease Mix	0.5 ml	80-6501-42
Protein G HP SpinTrap	16 columns	28-9031-34
Protein G HP MultiTrap	4 × 96-well plates	28-9031-35
NHS HP SpinTrap	5 ml medium and 24 empty spin columns	28-9031-28
Streptavidin HP SpinTrap	16 columns	28-9031-30
Streptavidin HP MultiTrap	4 × 96-well plates	28-9031-31
Ab SpinTrap	50 × 100 µl	28-4083-47
Ab Buffer Kit	1	28-9030-59

Literature

Antibody purification Handbook	18-1037-46
Affinity Chromatography Handbook	18-1022-29
Affinity Chromatography Columns and Media Selection Guide	18-1121-86

CyDye, Cy3, Cy5, Deep Purple, Drop design, Ettan, HiTrap, MultiTrap, Sepharose, and SpinTrap are trademarks of GE Healthcare companies. GE, imagination at work, and GE monogram are trademarks of General Electric Company.

CyDye or portions thereof is manufactured under license from Carnegie Mellon University under patent number 5268486 and other patents pending.

Deep Purple Total Protein Stain is exclusively licensed to GE Healthcare from Fluorotechnics Pty Ltd. Deep Purple Total Protein Stain may only be used for applications in life science research. Deep Purple is covered under a granted patent in New Zealand entitled "Fluorescent Compounds", patent number 522291 and equivalent patents and patent applications in other countries.

All goods and services are sold subject to the terms and conditions of sale of the company within GE Healthcare that supplies them. General Electric Company reserves the right, subject to any regulatory and contractual approval, if required, to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your local GE Healthcare representative for the most current information.

© 2006-2007 General Electric Company – All rights reserved.
 First published Nov. 2006.

GE Healthcare Bio-Sciences AB, a General Electric Company.

GE Healthcare Europe GmbH
 Munzinger Strasse 5, D-79111 Freiburg, Germany

GE Healthcare UK Ltd
 Amersham Place, Little Chalfont, Buckinghamshire, HP7 9NA, UK

GE Healthcare Bio-Sciences Corp
 800 Centennial Avenue, P.O. Box 1327, Piscataway, NJ 08855-1327, USA

GE Healthcare Bio-Sciences KK
 Sanken Bldg. 3-25-1, Hyakunincho, Shinjuku-ku, Tokyo 169-0073, Japan